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Szanowny kliencie

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NINGBO SEAGO ELECTRIC CO.,LTD

CUSTOMER SERVICE NUMBER:4008-843-006
ADDRESS:NO.239 Jianghu Road,Jiangbei District, Ningbo
WEB:www.seagoln.cn
EMAIL:service@seago.cn

SONIC TOOTHBRUSH
MODEL:S9


Dear Customer

Thank you for choosing seago sonic toothbrush. Seago is committed to provide high cost-effective products to you, and help you clean your teeth everyday to protect your oral health. A clinical trial report from a domestic well-known dental hospital shows,that continuously using of seago sonic toothbrush for 7 days,the plaque index will fell by 45%. It can help reduce the accumulation of plaque, smoke stains, tea and coffee stains, etc.


Important safeguards

- There are some people may not adapt such as stimulation and the feeling of pins and needles when using for the first time, please don't worry, this is the normal phenomenon, this phenomenon will fade out after several times using, and you will be comfortable to use this toothbrush in one week.
- Do NOT open, repair or alter the product without authorization.
- Do NOT knock or drop the product, if product is defective, please stop using.
- Please check if the local voltage is suitable with the marked voltage before connecting the power supply.
- If your oral or gum had surgery during the past two months, please consult your dentist before use.
- If you teeth is bleeding after using the toothbrush or bleeding continues more than one week, please consult your dentist.
- This product was designed to only clean your teeth, gums and tongue.
- Persons (including children) with a lack of experience and knowledge should not use this product, unless they were supervised to ensure proper usage.
- Do NOT use the charging base outdoors or high temperature. Do NOT put the base into the water or any liquid. Please cut off power supply when not in Use or before cleaning it. Please do not clean with hair care products (such as gel, oil), hand soap, or acid cleaning fluid otherwise it will damage the surface of the product.
- Children should be supervised to ensure proper usage and make sure they do not play with the toothbrush.
- Do NOT place in high temperature condition or directly sunning under the sunlight.

Component Name

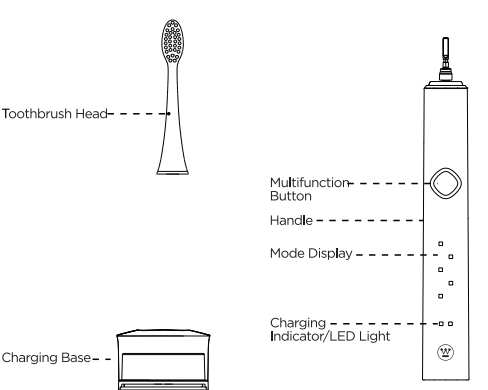Instrukcja

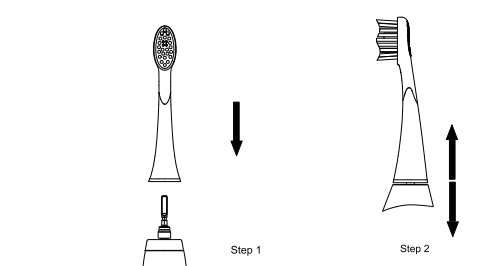Features

- With twin-engine magnetic levitation motor brings 32000-42000 micro-vibrations strokes per minutes, and it makes 6mm wide swing on the bristle to clean oral effectively.
- With 6 modes(White, Clean, Sensitive, Polish, Massage, Study) to meet different demands.
- The design of the fuselage accords with human body engineering. It's easy hold anti-slip and waterproof.
- Waterproof design and could meet the daily cleaning requirement enough.
- Safe and durable lithium battery. Full charged can last for 25 days, with power shortage reminder.
- Smart timer of 30s reminder and 2 minutes to educate people to brush teeth in the right way.
- High quality bristle of international brand promotes naturally clean teeth by gently polishing off surface stains. replaceable brush head.

How to use:

- Insert the toothbrush head tightly on the shaft lever.
- Wet the brush head with water and add toothpaste (the dosage is half less than ordinary toothpaste).
- Gently place the bristles of the brush on your teeth towards the gum line.
- Select the appropriate mode by pressing the mode selector button as shown in the diagram. The mode indicator button will show which mode is selected. Every 30 seconds will pause 0.5 seconds to prompt the user to change position in the process of brushing your teeth. The total time is 2 minutes. Toothbrush will stop working after 2 minutes. You can choose the different modes from the following.
 - a. **White** (32000/min Strong) - White mode can remove surface stains, to brighten your front teeth.
 - b. **Clean** (38000/min Middle) - Extra gentle mode for sensitive teeth and gums. Recommended for everyday use.
 - c. **Sensitive** (42000/min Weak) - Some people may feel uncomfortable when use at the first time. The bristles vibration strength changes from weak to strong to help you easily into the Sonic toothbrush experience.
 - d. **Polish** (32000/min Strong & Polish, 1 minutes) - Quick 2-minute mode to brighten and polish your teeth.
 - e. **Massage** (32000/min) - Massage mode to gently stimulate and massage your gums to improve gum health.
 - f. **Study** (32000/min Weak, the first 20 times increasing) - Someone don't fit it at the first time, who can begin with the mode.
- Gently move the brush head slowly across the teeth in a small back and forth motion to reach the total deeply cleaning effect. be sure to allow the bristles to reach between your teeth and massage your gums, please use the vibration strength from the toothbrush itself, do not force the bristles hard against the teeth during using.
- Clean the brush head and handle by washing and make sure toothbrush is cleaned.

How to charge:

- Inductive charging, users only need to put the toothbrush on the charging base which had connected power source. And toothbrush will be recharged.
- The indicator light will be on for about 21 seconds when charging and then turn off. It may takes about 16 hours for charging and the charging unit will turn off automatically when it is fully charged. And the indicator light will be on for about 21 seconds when the toothbrush is moved from the charging base.
- The multifunction button can't be used when charging.
- If I put handle in working on charging base, it will stop working and enter into charging mode.

Clean



- Handle
 - Clean the handle with damp soft cloth and air dry it after removing the toothbrush head.
- Do not keep the handle in water for long.
- Charging base
 - Clean the base with damp soft cloth and air dry it.
- Do not keep it in water.
- Toothbrush head
 - Clean and air dry it after use.
- Do not press or pull the head hard as it is easy to be deformed.
- Ensure cleanness effect of toothbrush, please change the toothbrush head every 3 months.

Trouble Shooting

Case	Possible cause	Solution
Vibration strength become weaken	1. Force the bristles hard against the teeth	1. Gently put the brush head on your teeth
The work time is short	1. The handle is not in the right place which causing the charging failure. 2. The battery maybe charge fail	1. Ensure to put it in a right way when charging 2. To extend the charging time
The product doesn't work	1. Checking if the product is full charged	1. Please contact the Seago customer service to deal with it.

Take out the battery

- This process is not reversible.
- Follow basis safety precautions, and protect your eyes, hands, fingers, and the surface of your work area.

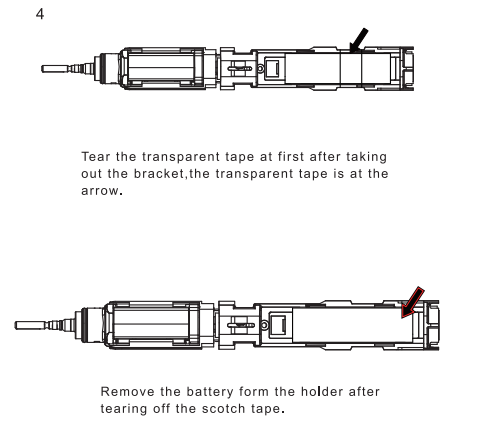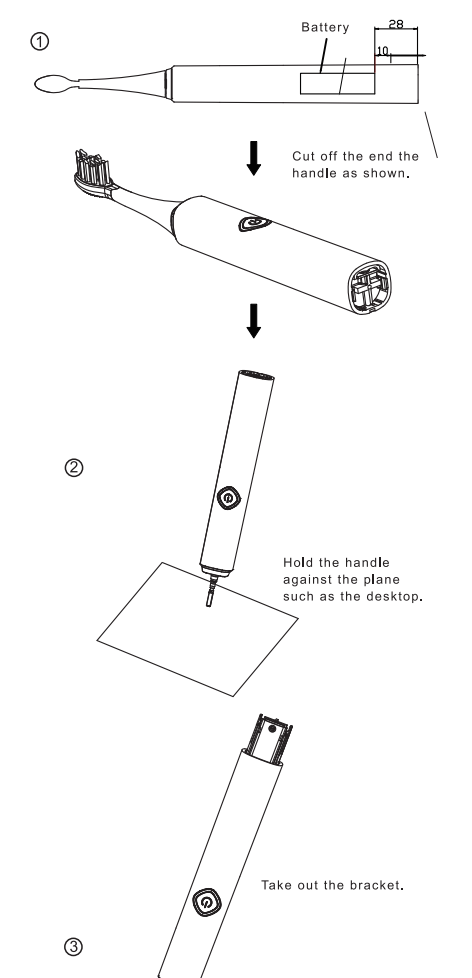Product information:

Product name: Sonic toothbrush	Model: S9
Input voltage: 5V	Input power: 1.5W
Working temperature: 0-50 °C	Battery capacity: 700mAh
Charging time: about 16 hours	Origin: China
Operative standard: GB4706.59-2008	

Environmental protection

Please don't discard the product that the service life is ending along with general household waste, please send it in the recycle bin.

Warranty Card

Thank you for choosing Seago sonic toothbrush and toothbrush sanitizers. All products have successfully passed strict quality assurance processes, Seago warrants the products to be free from defects under normal use and service for the original purchaser. The warranty period starts from the date printed on valid invoice for a period of one year, unauthorized disassembly or repair will void the warranty. The product can be returned within 7 days or exchange within 15 days after purchase.

Please fill below information by customer

Customer name _____ Purchasing date _____

Contact information _____

Please fill below information by retailer

Model number _____ Manufacture date _____

Retailer name _____

NINGBO SEAGO ELECTRIC CO.,LTD

CUSTOMER SERVICE NUMBER:4008-843-006
ADDRESS:NO.239 Jianghu Road,Jiangbei District, Ningbo
WEB:www.seagoln.cn
EMAIL:service@seago.cn